RO

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Warranty terms

Legal provisions on consumer protection are applicable.

1. Warranty period for this product:

1.1. 2 years for individuals.

1.2. 1 year for legal entities.

2. The warranty is calculated as of the date of delivery of the product.

During the warranty period, we will be responsible for repairing free of

charge any defects caused by the quality of the product itself under

normal conditions of use.

3. During the warranty period, if the defect is caused by the quality of the

product, please contact us.

4. Please keep your purchase invoice, as this serves as proof of the product

warranty.

5. The product loses its warranty in the following cases:

5.1. Failure caused by not using the appliance in accordance with the

instructions provided.

5.2. Damage caused by repairs by unauthorized persons.

5.3. Damage, scratches or breakage caused by moving or falling.

5.4. Damage caused by improper storage, maintenance or use by the user.

5.5. Damage caused by natural disasters or accidents.

5.6. Damage caused through the fault of the buyer (accidents including

transport)

5.7. Damage caused by falling, crushing, impact, exposure to water.

Do not remove screws and do not try to disassemble the appliance, otherwise the product warranty becomes null!

Recycling

Don't throw the device away!

Please contact a local organization to recycle this product at the

end of its useful life.

Product certification

This product complies with European conformity standards.

Technical specifications

Product Code	DET3DMAX	Message modes	3	Name	Shiatsu Back Massager with heat	
					Voltage	12V === 2.5A
Power	30 W	Timer	15 minutes			
	Voltage supply		100-240 V ~ 50/60 Hz	Dimension	114 x 46 cm	
Weight		Approx 5,2 kg				

Error	Reason	Method of resolution
<p>Device fails to start</p>	<p>Check if the device has been used for more than 30 minutes.</p> <p>1. If yes, allow the device to cool for 15 minutes before proceeding.</p> <p>2. The device may take longer to cool down when the ambient temperature is high.</p>	<p>Please contact your retailer or an authorised service centre if none of the above methods solve the problem.</p>
<p>Device fails to start</p>	<p>Device in not ON</p>	<p>Press the Power ON button (1)</p>
<p>Device suddenly stops</p>	<p>The automatic switch-off function is activated after 15 minutes of operation.</p>	<p>1. Continue using by pressing the "Power/ ON" button. After every 30 min of use switch off the appliance for 15 minutes.</p>
<p>There is a clattering noise during use</p>	<p>This is normally caused by the operation of the motor and rollers during use.</p>	

- Before cleaning, make sure the appliance is switched off and the adapter has been unplugged.
- Avoid immersing the appliance in any liquid during cleaning.
- Clean only with a soft, slightly damp sponge.
- Do not allow water or any other liquid to get inside the unit.
- Avoid using abrasive cleaners, brushes, gasoline, kerosene, glass/ furniture polishes.

Cleaning instructions

If you do not use the appliance for a long period of time, please keep it in a safe, dry and covered place to prevent dust accumulation, out of direct sunlight and out of reach of children.


Avoid contact with sharp edges or sharp objects that could cut or puncture the surface.

Store the power adapter, after coiling the cord, together with the appliance. Do not hang the appliance by the power cord.

Maintenance and storage

Instructions for use

Use the remote control below to customize your back massage experience:

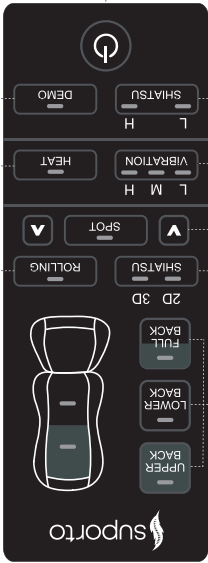Massage Area Buttons: Select upper, lower or full back to designate the area covered by the massage.

Shatsu Button: Press to shift 2D or 3D shatsu style.

Spot Button: Press to keep Shatsu function in one area, Use the arrow buttons to adjust the position up or down.

Vibration Button: Press to activate the seat vibration. You may cycle through low, medium, high and off settings.

Shatsu Intensity Button: Press to set the shatsu intensity at the low or high level.

Power Button: Press to turn the device on or off.

Rolling Button: press to activate the action rolling to travel up and down across the muscle gently. Continue pressing the button to adjust the width of the massage rollers to custom fit the body.

Heat Button
Press to activate heat on shatsu nodes. The heat can not be turned on if the massage function is not activated.

Demo Button: Quick demonstration of the lights will blink to indicate the active functions. The L.E.D. device will automatically turn off. Note : While Demo function be activated, the other functions can not be turned on.

day.


You should allow your body some rest periods. To avoid straining your muscles too much, we recommend not exceeding 15 minutes of continuous massage per session and not more than 3 sessions per day.

Attention !

The device will automatically switch off after 15 minutes (timer run out).

If the Shatsu can't reach your lower back area: If you want to massage your lower back, but the shatsu nodes can not reach that low. You may adjust the height of your back when sitting, by placing a seat cushion or a pillow on the sitting part.

If the Shatsu is too strong: Install the flap onto the cushion to get a softer massage. Place a blanket or towel between you and the massager to achieve a far softer massage intensity if you still feel too strong shatsu after installing the flap onto the cushion.

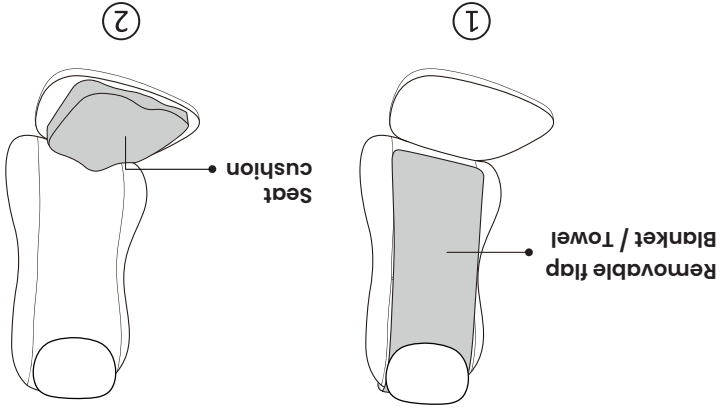Correct use of the device

Make sure that the device and accessories are free from visible defects

before use.

If in doubt, do not use the device and contact the retailer or an authorized

service center.


Follow the instructions described below only after you have read, understood and observed all the safety rules and measures described in this document.

Setup

1. Attach the massager to a chair with elastic straps or place it on other

support.

2. Connect the adapter cable to the corresponding cable on the cushion.

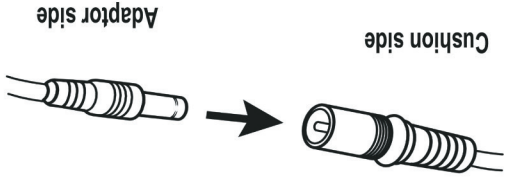3. Plug the home adapter into an electric outlet.

4. Turn on the device using the controller (instructions on page 8).

5. When finished, place the controller in the pouch at the side of the cushion.

6. You also may soften the massage pressure by placing a blanket or towel between you and the device.

3. During sleep. Be careful not to fall asleep while the machine is still working.
 2. Near flammable or corrosive substances.
 1. Near fire or a heat source.
- Do not use in the following cases:

Fire risk instructions

8. Do not pull, twist or bend the power cord. Avoid placing the power cord over or in contact with sharp or pointed objects.
 7. Do not repair the appliance yourself, contact the retailer or an authorized service center. Improper repairs may expose the user to serious risk of injury and the product warranty can become null.
 6. If the appliance has been dropped or otherwise damaged as a result of being dropped, it may be used only after being checked by authorized personnel.
 5. If the appliance has been dropped or otherwise damaged as a result of minutes after every hour of use.
 4. Do not use the appliance for more than 30 minutes consecutively, there is a danger of overheating. Allow the appliance to cool for at least 30 minutes after every hour of use.
 3. Handle only with dry hands. Do not unplug the power plug with wet hands.
 2. Keep the power cords away from hot surfaces.
 1. Do not pull the power cord or the appliance to unplug the plug. Never hold the appliance by the power cord.
 - 1.5. If water has entered the appliance, the power cord, power plug or anywhere inside the product.
 - 1.4. In the presence of sources of increased humidity (bath, sauna).
 - 1.3. During storms.
 - 1.2. (burningsmell, broken or damaged wire, visible burns, physical defects, etc.).
 - 1.1. At a mains voltage different from that stated on the appliance.
1. Do NOT use in the following cases:

Only use the power adapter supplied with the device.

Since it is an electrical device, it shall be used with care and caution to avoid the risk of electric shock.

The appliance is very well insulated against the risk of electric shock.

Measures regarding the risk of electric shock

9. Stop using the product if you notice malfunction or any other abnormal operation.
10. Do not allow the massage cushion to come into contact with pointed or sharp objects.

Ask your doctor advice before use

- Ask your doctor advice before use if you are in at least one of the following cases:
1. People with chronic back conditions, who have recently undergone surgery, open wounds, swelling, trauma, burns to the body etc.
 2. People with diabetes.
 3. People with pain whose cause is unknown.
 4. People with psychiatric conditions.
 5. People with pacemakers or implanted artificial organs.
 6. People with heart dysfunction, tumors and hypertension.
 7. People with brain haemorrhage and thrombosis.
 8. People with haemorrhagic diseases, septiccaemia, fever, etc.
 9. Women who are pregnant or menstruating.
 10. People with advanced osteoporosis or spinal fractures.
 11. People with heat intolerance.
- Please also consult your doctor** in other cases if you are not sure whether the device is right for you or if you are undergoing treatment prescribed by your doctor.
- Avoid using the product before bedtime, it can affect your sleep.
- If you experience pain or discomfort during the session, stop using the device and consult your doctor!

Safety rules and measures

General safety measures

1. This appliance is NOT a toy. Close supervision is required when this appliance is used near children or by persons with disabilities. Do not allow children to play with the appliance.
2. The massage device must only be used for its intended purpose as described in the instruction manual.
3. To avoid excessive stimulation of muscles and nerves, the recommended duration of massage should not exceed 15 minutes per session. Allow a waiting time between sessions.
4. Do not use the machine in damp or dusty environments as this may lead to malfunction or electric shock.
5. The appliance is not intended for commercial use, but for individual use only.
6. Do not sit on the appliance during use.
7. Protect from excessive heat or direct sunlight to prevent damage.
8. Do not use the appliance if the fabric cover is damaged.

Importance of reading instructions

Taking the time to read these instructions is an investment in your safety and well-being.

The back massager is an electrical product! Read these operating instructions carefully!


Failure to observe the following instructions may result in property damage and personal injury. Keep the operating instructions and present them to other users. Hand in this manual when you dispose of the appliance.

For instructions in online version, go to this link or scan the QR code:


<https://www.supor.ro/blog/instructiuni/masaj-spate>


Package content

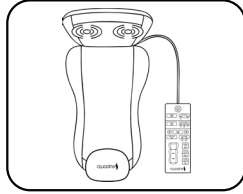1 x Back massager
+ remote control

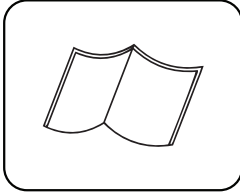1 x User manual

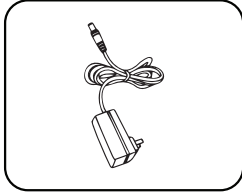1 x Power cord adapter

This back massager features 10 flexible 2D/3D nodes combined with a vibrating cushion to deliver deep relaxation and effective muscle relief. With 2 intensity levels and customizable options—Spot, lumbar, or upper back—you can tailor the massage exactly to your needs.

Soothing infrared heat boosts circulation and eases stiffness, while two massage modes (Demo and Roll) add comfort and variety. The ergonomic design, with a soft cushion cover and secure elastic chair strap, ensures stability and convenience. Plus, the built-in 15-minute auto timer makes every session safe and worry-free.

Device presentation

“Take care of your body,
it's the only place you have to live.”

Jim Rohn

Thank you for buying this product!
We value your trust in Supporto.
Relax and improve your health with this
massage device!

Thank you	3
Device presentation.....	3
Importance of reading instructions.....	4
Package content.....	4
Ask your doctor advice before use	5
Safety rules and measures.....	5
General safety measures	5
Measures regarding the risk of electric shock.....	6
Fire risk instructions	6
Correct use of the device	7
Setup.....	7
Operating instructions.....	8
Maintenance and storage.....	10
Troubleshoot.....	11
Technical specifications.....	12
Warranty terms	13
Recycling.....	13
Product certification.....	13

Detenso 3D Max - Shiatsu Back Massager with heat

USERS MANUAL

This product complies with European
conformity standards.

Made to order
Supporto Ortopedic SRL
djtutor@supporto.ro
+40757246463
Importer: Supporto Ortopedic SRL
www.supporto.ro