RO

Pistol de Masaj F5 Titan

MANUAL DE UTILIZARE

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Jim Rohn

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The Massage Gun is a portable device that uses percussion therapy to soothe muscles and enhance recovery.

Equipped with a titanium alloy massage head, a powerful motor, and multiple speed settings, the Suporto F5 Titan is the perfect choice for athletes and fitness enthusiasts seeking a professional-grade massage.

With its ability to target deep muscle layers, it effectively prepares muscles before exercise and helps relax and relieve muscle tension afterward.

Its long-lasting battery and ergonomic design make it a dependable companion for consistent muscle care.

Whether you're dealing with fatigue, aiming to prevent discomfort, or in need of recovery, the F5 Titan Support is an essential tool for your health and well-being.

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Warranty terms

Legal provisions on consumer protection are applicable.

Do not remove screws and do not try to disassemble the appliance, otherwise the product warranty becomes null!

1. Warranty period for this product:

1.1. For individuals 2 years.

1.2. For legal entities 1 year.

2. The warranty is calculated as of the date of delivery of the product.

During the warranty period, we will be responsible for repairing free of charge any defects caused by the quality of the product itself under normal conditions of use.

3. During the warranty period, if the defect is caused by the quality of the

product, please contact us.

4. Please keep your purchase invoice, as this serves as proof of the product

warranty.

5. The product loses its warranty in the following cases:

5.1. Failure caused by not using the appliance in accordance with the

instructions provided.

5.2. Damage caused by repairs by unauthorized persons.

5.3. Damage, scratches or breakage caused by moving or falling.

5.4. Damage caused by improper storage, maintenance or use by the user.

5.5. Damage caused by natural disasters or accidents.

5.6. Damage caused through the fault of the buyer (accidents including

transport)

5.7. Damage caused by falling, crushing, impact, exposure to water.

Product certification

This product complies with European conformity standards.

Technical Specifications

Product type	Massage gun
Model	FS Titan
Voltage	7.3V~
Amperage	2A
Power	25W
Amplitude	10 mm
Speed rpm	2000-3200
Product size	13.6 x 6.1 x 19.3 cm
Weight	0.680 kg
Charging time	~ 180 min
Battery capacity	4000 mAh
Timing	10 minutes
Operating mode	ON/OFF Button
Massage gears	5

Transport Conditions:

- The product must be stored in the original box upright, carefully, lightly and protected from rain and moisture.
- Ambient temperature: 20 °C ~ +35 °C
- Relative humidity: ≤ 90%

Environmental conditions for product

Installation:

- Ambient temperature: +5 °C ~ +40 °C
- Relative humidity: ≤ 90%
- Atmospheric pressure: 70 kPa ~ 106 kPa

Recycling

falls:

- This product is suitable for transportation in automobiles, trains, ships, airplanes and other modes of transportation.
- During transportation, the product should be secured to prevent violent vibration, collisions, overturning and falls;
- Do not throw the product in the trash with household waste!
- To prevent possible harm to the environment or human health through uncontrolled waste disposal, recycle this product responsibly helping to reduce pollution and conserve natural resources.
- Please contact a local organization to recycle this product at the end of its useful life.
- The product will be recycled after disposal of the batteries.

Error management

- Battery removal procedure (for authorized personnel):**
1. Ask the seller for disassembly instructions for this appliance.
 2. Disconnect the appliance from the charger.
 3. Use professional tools to open the plastic housing.
 4. Cut the battery cables with scissors.
 5. Remove the battery.

Problem	Cause	Solution
The device won't start	The power is off. The battery is discharged.	Start the device by pressing the On/Off button firmly. Charge the device.
The device turns off automatically	Overheat protection is activated. The appliance will automatically switch off when overloaded. The device is set to turn off automatically after 10 minutes.	Restart the device after a 10-minute break. Restart the device for another 10 minutes of use.
The device is locked	The lubricant solidifies due to the low temperature.	Gently loosen the message head mounting slot, then turn the device on.
The message head falls out	The message head is loose or the rubber gasket is worn.	Install the message head as instructed or change the rubber gasket.

Battery disposal and recycling

1. At the end of the operating period the battery must be safely disposed of and recycled. Do not dispose of the battery in household waste!
2. Please support the recycling of rechargeable batteries. The lithium-ion batteries used in this product are a valuable resource that can be recycled, so please **dispose** of them at a **specialized service** and then hand them over to a **professional organization that can recycle rechargeable batteries**.
3. The battery disposal procedure will only be performed by **authorized personnel**. Failure to follow these instructions may result in fire, electric shock, injury and other accidents!
4. Consult the seller for any repair requests.

Battery precautions

1. Do not disassemble the battery in any way (except authorized pers.).
2. Keep out of the reach of infants and small children.
3. The battery in this product is a special rechargeable battery and should not be used in other products.
4. Battery life may vary depending on conditions of use, charging and storage.
5. Keep the battery away from water, fire, and heat sources (e.g. radiators, stoves, fireplaces, fireplaces, stoves, direct sunlight, etc.).
6. Do not put the battery in microwave ovens or pressure cookers.
7. Do not pierce the battery with sharp objects (nails, blades, screwdrivers, etc.).
8. Do not hit, slam, or stomp the battery.
9. Do not connect positive and negative battery terminals directly to metal.
10. Do not store the battery in the same place where you keep necklaces, hair clips.
11. Do not use the battery if you notice any abnormality such as leakage, odor, heat, deformation, discoloration or any other unusual condition.
- If the device is in use or charging, stop use immediately and disconnect the charger.
12. Avoid contact of the liquid drained from the battery with the skin, eyes or other parts of the body, it causes irritation and injury. In case of contact, wash with plenty of water and seek medical advice.

- **Clean** the device with a slightly damp cloth and wipe it with a soft cloth or dry wipes.
- Avoid using alcohol, thinner, benzene and other solvents. These can cause the device to malfunction, damage the casing and discoloration.
- When **storing** or travelling, place it in the carrying case, away from liquids, dust, frost or heat sources.
- Do not wrap the charging cord around the device when storing to avoid damage.
- These massage guns should be stored in a well-ventilated, dry place, free from corrosive gases, at a height of 25 cm above the floor.

Maintenance and storage

*If LED1 flashes the battery is below 10% and needs charging

battery charge level as follows:

- The device can operate approximately 2 to 4 hours after being fully charged, depending on the speed and pressure applied during use.
- During use, the LED indicator at the base of the handle indicates the battery charge level as follows:
 - It is not recommended to allow the battery to fully discharge.
 - Charge the device at least once a month if not in use to prolong battery life.
 - The device can operate approximately 2 to 4 hours after being fully charged, depending on the speed and pressure applied during use.
 - Do not leave the product charging overnight or unattended during charging or use.
 - It is not recommended to allow the battery to fully discharge.
 - Charge the device at least once a month if not in use to prolong battery life.
- After the first charge, the battery can be charged at any time regardless of the battery level.
- Do not leave the product charging overnight or unattended during charging or use.
- It is not recommended to allow the battery to fully discharge.
- Charge the device at least once a month if not in use to prolong battery life.
- The device can operate approximately 2 to 4 hours after being fully charged, depending on the speed and pressure applied during use.
- During use, the LED indicator at the base of the handle indicates the battery charge level as follows:

Device charging

Failure to follow the following charging instructions will result in premature battery wear and may cause property damage and injury. We are not liable for damages caused by failure to follow these instructions.

5V ⚡ 2A

Use with 5V ⚡ 2A
charger output

1. To charge the product, connect it to a **5V-1A or 5V-2A** power supply using the USB-A to USB-C cable supplied with the product.

2. Do not use the charging cable if it is damaged or poorly connected. It may cause electric shock, short circuit and fire.

3. The device can be charged by: universal charger, phone charger, external battery or laptop/computer (slower charging).

Charging time will be longer if the power supply is less than 5V-2A.

Do not charge the device with very powerful chargers at a higher voltage than indicated to avoid damage to the battery and risk of fire, explosion,

property damage and personal injury!

4. Charge the device fully before first use.

As long as the battery indicator lights light up in succession, the device is

charging.

The charging time varies depending on the power supply and battery charge

level, between 1 and 3 hours.

When the battery indicator lights stop blinking and remain fully lit, the device

is fully charged.

Do not turn the device on while charging!

5. To get the best battery capacity quickly, repeat a full charge – discharge

cycle 3 times.

6. Disconnect the charging cable by grasping the plastic/mold front end and


do not pull on the cable. In addition, do not disconnect with wet hands to

avoid electric shock, short circuit or fire.

1. Press and hold the On/Off button  on the panel for 3 seconds to start the unit in 1st Gear.
2. Press the button briefly to change between Gears 0-5:
3. When operating on Gears 1-5, press and hold the button to move to Gear 0, at which time the appliance will turn off (standby);
4. To turn the appliance completely off, return to Gear 0 then press and hold the button for 3 seconds.
- Left on Gear 0 (standby) for 20 seconds, the machine will automatically turn off if no operation is performed.
- After every 10 min. of use, the machine will automatically switch to Gear 0 and will stop if a new massage session is not started.
- The LED indicator above indicates the speed/power gear, Gear 1 being the lowest and Gear 5 being the highest.
- To massage effectively choose the right massage head and apply appropriate speed and pressure to each area of the body.

Correct use of device

2. To replace a massage head remove the massage head from the installation slot by pulling gently to detach. Then install the new head following the instructions in step 1.
- Caution!** The area where the massage head is installed may become hot due to continuous movement, please be aware of burn hazards when replacing the massage head and allow the unit to cool down if it is hot.

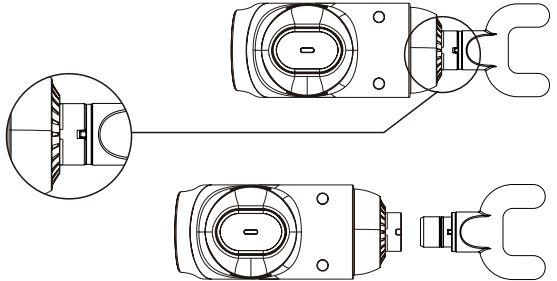Installing and replacing massage heads

Before installing or replacing a massage head, make sure the device is turned off.

1. To install a massage head, insert the head stem into the slot then push with light pressure until the rubber gasket meets the slot.

Warning! The flat titanium head screws on.

Make sure you have secured the head correctly!

Some massage heads have a protrusion in the rubber gasket area that

prevents the head from rotating during massage. This should fit snugly into

the slot as shown in the following picture.

Ball head: For general use, large muscle groups (chest, back, legs, quadriceps, shins, gluteals), small muscle groups (biceps, triceps, shoulders), for sensitive areas (cervical, neck and shoulder area, forearm, abdomen, paravertebral muscles) as well as for relieving muscle and joint pain.

Head with air cushion: For relaxing the muscles. General use, back, but especially for sensitive areas of the body, ankles, wrists.

U-shaped head: For back, neck, shoulders, muscles and soft, sensitive tissues and large muscle groups. **Warning: Do not touch the spine!** Back massage under the guidance of a specialist, otherwise there is a risk of serious injury.

Bullet-shaped head: For joints, deep tissues, small areas of muscles such as ankles and wrists, paravertebral muscles, trigger points, knots and muscle contractures.

Caution: Sensitive areas should be massaged under the guidance of a specialist.


Titanium flat head: For the whole body, muscle relaxation, except sensitive areas. **Avoid bony areas and areas near the spine!**

Electrical and mechanical precautions

1. Do not use the device for more than 30 min continuously to avoid overheating and premature wear. Allow a 10 minute break every 30 min of continuous use. Stop use immediately if it has overheated or emits a burning odor.
2. Do not personally open or service the appliance. Doing so exposes you to danger and voids the product warranty.
3. Do not use the appliance in the bath or shower. Do not put it in water, keep it away from liquids and areas of high humidity, and do not allow water to enter the vents and interior.
4. Do not drop, hit or knock the appliance. Do not use it if it malfunctions or is damaged due to dropping.
5. Do not charge the device while it is operating.
6. Do not use, charge or store the device near fire, flammable substances and heat sources.
7. Do not use the product near long hair and jewelry (chains, bracelets, piercings, etc.), on loose clothing and clothing with wide sleeves. Risk of entrapment, breakage and injury.
8. Do not use the product in areas where aerosol (spray) products are used or where oxygen is used.

Product description and parts

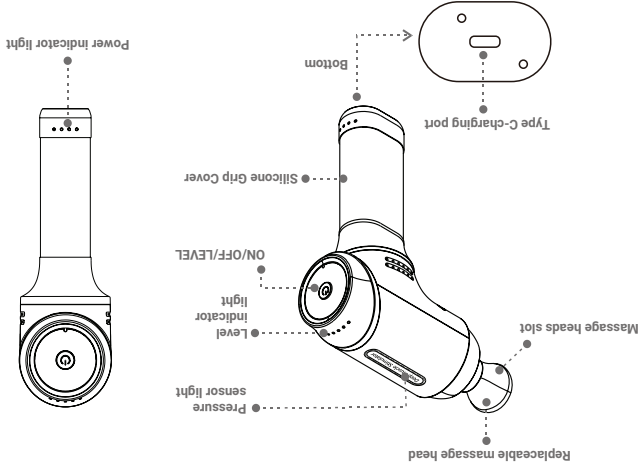Safety rules and measures

7. Check that the appliance starts and operates normally. If you notice abnormal operation, strange noises, smoky or burning smells, turn the appliance off immediately and contact the supplier.

8. A grease will be noticeable around the massage head installation slot. This is normal because the part is subject to friction and is heavily greased from the factory. Lightly wipe off excess grease with a dry tissue, but not completely, the area should remain greased.

1. The massage gun is intended for healthy adults only.
2. This product is not a toy. Keep out of the reach of children!
3. Seniors, children and people with limited cognitive abilities or sensory impairments should only use this product under supervision!
4. Use this product with care: Applying excessive force or speed during use may cause bruising.
5. Use only on soft tissue areas and avoid hard or bony areas of the body, including the head, cervical spine, spine, joints, shins, dorsal side of the legs and hands, etc.
6. Use only on dry, clean skin. Apply gentle pressure massaging about 60 seconds per muscle.
7. Check the massaged area frequently and stop immediately at the first sign of pain or discomfort. If pain or discomfort persists and cannot be relieved, seek medical advice immediately.
8. It is recommended to start the massage at the lowest intensity. Finish each session in 10 minutes. Avoid massaging in the same spot for more than 3 minutes, muscles and nerves can over-stimulate and you may injure yourself.
9. Do not use around sensitive parts, organs and central nervous system such as head, neck, spine, eyes, teeth, reproductive organs, artificial limbs, etc.
10. Do not use the massage gun on the genital area (e.g. for sexual purposes). This is a powerful device and can cause serious injury.
11. Do not apply the product to any area that may be fractured or broken or to areas where metal bone rods are installed.
12. Do not use the product on and/or near areas with implants (breast implants or prostheses) as it may weaken the implant or injure you.
13. Keep fingers, hair or other body parts away from the shaft and back of the vibrator head as there is a risk of pinching them.
14. Do not use on animals.

To reduce the risk of electric shock, fire and personal injury, this product should be used in accordance with the following measures:

Ask your doctor advice before use

DO NOT use the device without consulting your doctor if you are in any of the situations below:

1. People who have varicose veins, thrombotic diseases, dermatitis of any kind, skin infections, contusions, local lesions and bleeding, open wounds on the skin surface, venous stasis and rashes.
2. Persons with psychiatric conditions.
3. People suffering from diabetes and complications such as neuropathy or retinal disorders.
4. Women during menstruation or during pregnancy.
5. People suffering from heart disease, heart failure, using pacemakers or defibrillators, or those who have had recent surgery.
6. Persons suffering from malignant tumors, epilepsy, migraine headaches, cachexia, septicemia, herniated disc, spondylolisthesis, osteoporosis, osteophytes, fractures, influenza, fever.
7. Persons who have had recent joint replacement surgery or placement of intraterrine device, metal implants.
8. Persons advised by their doctor not to massage.
9. Other health problems unless you know for certain that the appliance is suitable.

Before first use


1. Open the box and check that the main device and accessories are in good condition. If you notice any damage, signs of deterioration or missing parts, contact the supplier immediately.
2. Do not leave packaging material, i.e. plastic bags, within the reach of children. There is a risk of suffocation.
3. Make sure the battery is fully charged. (see "Charging" at page 10)
4. Do not use the device while it is charging!
5. Make sure the ventilation holes of the appliance are not blocked.
6. Pick the desired massage head and connect it to the device. (See page 9)

Importance of reading instructions

Taking the time to read these instructions is an investment in your safety and well-being. The massage gun is an electrical product! Read these operating instructions carefully!

Failure to observe the following instructions may result in property damage and personal injury. Keep the operating instructions and present them to other users. Hand in this manual when you dispose of the appliance. For instructions in online version, go to this link or scan the QR code:

INSTRUCTIONS: <https://www.supporto.ro/blog/instructioni/massage-gun>

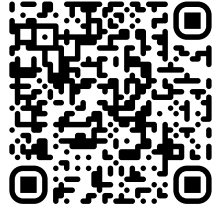Package Content

1 x Massage Gun

1 x User manual

1 x Cable power supply Type C

5 x Massage heads

1 x Massage head rubber gasket

The Massage Gun is a portable device that uses percussion therapy to relax muscles and improve muscle recovery.

Equipped with a titanium alloy massage head, powerful motor and multiple speed settings, the Supporto FS Titan is the ideal choice for athletes and fitness enthusiasts looking for a professional massage.

Able to penetrate into deep tissues, it is effective in preparing muscles before training and after by relaxing and relieving muscle tension. Its long-lasting battery life and ergonomic design make it a reliable partner for long-term muscle care.

Whether you're experiencing fatigue, want to prevent discomfort or need to recover, the FS Titan Supporto is a necessary tool for your health and well-being.

Device presentation

“Take care of your body,
it's the only place you have to live.”
Jim Rohn

Thank you for buying this product!
We value your trust in Supporto.
Relax and improve your health with this
massage device!

Contents

3	Thank you
3	Device presentation
4	Importance of reading instructions
4	Package Content
5	Ask your doctor advice before use
5	Before first use
6	Safety rules and measures
7	Electrical and mechanical precautions
7	Product description and parts
8	Installing and replacing massage heads
9	Correct use of the device
10	Device charging
11	Maintenance and storage
12	Battery disposal and recycling
12	Battery precautions
13	Error Management
14	Technical specifications
14	Recycling
15	Warranty terms
15	Product certification

Made to order
Support Ortopedic SRL
Support Ortopedic SRL
djutor@suporto.ro
+40757246463
Importer: Support Ortopedic SRL
www.suporto.ro

This product complies with European
conformity standards.

USER MANUAL

Massage Gun FS Titan

EN