RO

APARAT MASAJ PICIOARE UV

MANUAL DE UTILIZARE

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Warranty terms

Legal provisions on consumer protection are applicable.

1. Warranty period for this product:

1.1. 2 years for individuals.

1.2. 1 year for legal entities.

2. The warranty is calculated as of the date of delivery of the product. During the warranty period, we will be responsible for repairing free of charge any defects caused by the quality of the product itself under normal conditions of use.

3. During the warranty period, if the defect is caused by the quality of the

product, please contact us.

4. Please keep your purchase invoice, as this serves as proof of the product warranty.

5. The product loses its warranty in the following cases:

5.1. Failure caused by not using the appliance in accordance with the instructions provided.

5.2. Damage caused by repairs by unauthorized persons.

5.3. Damage, scratches or breakage caused by moving or falling.

5.4. Damage caused by improper storage, maintenance or use by the user.

5.5. Damage caused by natural disasters or accidents.

5.6. Damage caused through the fault of the buyer (accidents including transport)

5.7. Damage caused by falling, crushing, impact, exposure to water.

Do not remove screws and do not try to disassemble the appliance, otherwise the product warranty becomes null!

Recycling

Don't throw the device away!
Please contact a local organization to recycle this product at the

end of its useful life.

Product certification

This product complies with European conformity standards.

Technical specifications

Product Code	AMP2N	Message modes	4
Name	Foot Massager UV	Mode of operation	Buttons
Voltage	12V~	Display mode	Illuminated panel
Power	36W	Heating mode	Infrared heating with carbon fiber pads
Amperage	3A	Sterilization method	UV Light
Operating temperature	0-40°C	Power supply	Adapter: Input 100-240V AC, 50/60Hz, Output 12V DC
Weight	4.2 kg	Dimensions	320 x 356 x 210 mm

Error	Reason	Method of resolution
Device fails to start	Not connected to the socket	1. Plug the adapter into the socket. 2. Check if the adapter plug is connected to the charging port of the device. 3. Check that the socket has power
	Device in not ON	Press the Power ON button (1)
	Check if the device has been used for more than 30 minutes.	1. If yes, allow the device to cool for 15 minutes before proceeding. 2. The device may take longer to cool down when the ambient temperature is high. Please contact your retailer or an authorised service centre if none of the above methods solve the problem.
Device suddenly stops	The automatic switch-off function is activated after 15 minutes of operation.	1. Continue using by pressing the "Power/ ON" button. 2. After every 30 min of use switch off the appliance for 15 minutes.
	There is a clattering noise during use	This is normally caused by the operation of the motor and rollers during use.

Maintenance and storage

Before cleaning, make sure the machine is switched off and the adapter has been removed from the socket.
If you do not use the appliance for a long period of time, please keep it in a safe, dry and covered place to prevent dust build-up, out of direct sunlight and out of the reach of children.
The appliance can be wiped clean with a slightly damp cloth on the outside, without detergent.

Cleaning instructions for removable covers

The covers are removed by unzipping the zipper and unfastening the Velcro inside.

Hand wash at 30 degrees C or automatic wash at 30 degrees C but without spinning.


No bleach

No ironing

Washing
Temperature

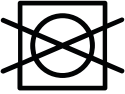No spin drying

Covers material: 95% polyester; 5% spandex.

Operating Instructions

Please choose a chair of the right height and sit down, then put your feet into the massager. Your feet should be relaxed, in a comfortable, unstretched position.

Do not stand in the device, it will get damaged!

Do not set up and use the device in a bathroom or similar/wet area.

Operate according to the control panel below to achieve the desired foot massage:

1. **On/Off button:** Press to start the machine. It will start in standard Mode C
- Relaxation massage for 15 minutes. The active functions are illuminated blue on the panel. Press again to switch off the machine.
2. **Massage Modes key:** Press to switch between massage types A, B, C and D. The corresponding program indicator will light up:
 - Mode A: Tension reduction,
 - Mode B: Tired feet,
 - Mode C: Comfortable and relaxing,
 - Mode D: Vitality to improve sleep.

3. **Compression Key:** Press to start the compression function, press again to change the intensity level: low/high or off.

4. **Heat key:** Press to turn on the Heat function, press again to change heat level: low/high or off.


5. **Rolling/Kneading key:** Press to turn on the kneading massage function, press again for low/high or off.

6. **Timer key:** Press to switch the timer between 15min or 30min. The corresponding indicator will light up.

7. **UV Light key:** Ultraviolet sterilisation function activated. Sterilization

starts only after putting the top lid back on after the massage session. The "UV" indicator light on the lid remains lit until sterilization is complete.

Product description and parts

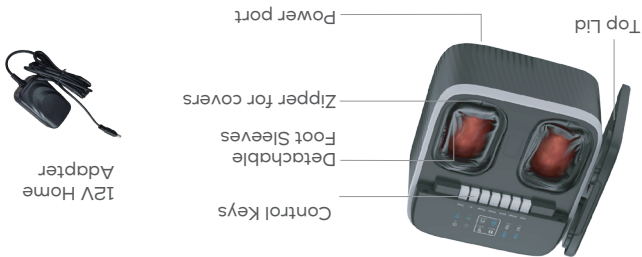*The picture is for reference only

Correct use of the device

Make sure that the device and accessories are free from visible defects before use.

If in doubt, do not use the device and contact the retailer or an authorized service center.

Follow the instructions described below only after you have read, understood and observed all the safety rules and measures described in this document.

Setup

1. Place the device where you want to use it. The surface must be flat and dry.
2. Connect the cable from the adapter to the power port of the appliance.
3. Plug the adapter into the electrical outlet (the adapter is safe to use at 110-120V & 220-240V)
4. After plugging into the outlet, execute the commands using the buttons as described below.

3. During sleep. Be careful not to fall asleep while the machine is still working.
 2. Near flammable or corrosive substances.
 1. Near fire or a heat source.
- Do not use in the following cases:

Fire risk instructions

9. Do not pull, twist or bend the power cord. Avoid placing the power cord over or in contact with sharp or pointed objects.
8. Do not repair the appliance yourself, contact the retailer or an authorized service center. Improper repairs may expose the user to serious risk of injury and the product warranty can become null.
7. If the appliance has been dropped or otherwise damaged as a result of being dropped, it may be used only after being checked by authorized personnel.
6. Do not use the appliance for more than one hour consecutively, there is a danger of overheating. Allow the appliance to cool for at least 30 minutes after every hour of use.
5. Handle only with dry hands. Do not unplug the power plug with wet hands.
4. Keep the power cords away from hot surfaces.
3. Do not pull the power cord or the appliance to unplug the plug. Never hold the appliance by the power cord.
2. If you notice a malfunction, switch off the appliance immediately and unplug it. Look for any damage to the product.
- 1.5. If water has entered the appliance, the power cord, power plug.
- 1.4. In the presence of sources of increased humidity (bath, sauna).
- 1.3. During storms.
- 1.2. If the appliance or its accessories have any noticeable defects (burning smell, broken or damaged wire, visible burns, physical defects, etc.).
- 1.1. At a mains voltage different from that stated on the appliance.

1. Do not use in the following cases:

the risk of electric shock:

Since it is an electrical device, it shall be used with care and caution to avoid the risk of electric shock.

Measures regarding the risk of electric shock

16. Do not sit or store heavy objects on the device.
15. Do not store the product in direct sunlight to avoid discoloration over time.
14. Pull carefully on the zipper when removing the foot cover to avoid damaging it.
13. Place the power cord in an area where there is no danger of tripping.
12. Do not stand while your feet are in the appliance to avoid damaging the appliance and causing danger.

Ask your doctor advice before use

- Ask your doctor advice before use if you are in at least one of the following cases:
1. People with chronic foot conditions, who have recently undergone foot surgery, open wounds, swelling, trauma, burns to the feet, etc.
 2. People with diabetes.
 3. People with pain whose cause is unknown.
 4. People with psychiatric conditions.
 5. People with pacemakers or implanted artificial organs.
 6. People with heart dysfunction, tumors and hypertension.
 7. People with brain haemorrhage and thrombosis.
 8. People with haemorrhagic diseases, septicæmia, fever, etc.
 9. Women who are pregnant or menstruating.
 10. People with advanced osteoporosis or spinal fractures.
 11. People with heat intolerance.
- Please also consult your doctor in other cases if you are not sure whether the device is right for you!
- For children and people with disabilities use is allowed under adult supervision!
- If you experience pain or discomfort during the session, stop using the device and consult your doctor!

Safety rules and measures

General safety measures

1. The appliance must be switched off and unplugged after each use and before each cleaning operation.
2. Do not leave the appliance unattended when children are nearby.
3. This appliance is NOT a toy. Careful supervision is required when the appliance is used by or near children or persons with disabilities. Do not allow children to play with the appliance.
4. To avoid over-stimulation of muscles and nerves, the recommended duration of massage should not exceed 30 minutes per session. Allow a waiting time between sessions.
5. Use the appliance only at ambient temperatures of 0-40°C.
6. The appliance is not intended for commercial use, for individual use only.
7. Do not put pressure on the upper area of the appliance or on buttons.
8. Do not leave packaging material within reach of children. There is a danger of asphyxiation.
9. Use only in enclosed space.
10. Do not use it on animals.
11. Do not use it for any purpose other than foot massage.

Importance of reading instructions

Taking the time to read these instructions is an investment in your safety and well-being.
The foot massager is an electrical product! Read these instructions for use carefully!


Failure to follow these instructions may result in property damage or personal injury. Keep the instructions for use and show them to other users. Hand this manual in when you dispose of the appliance.

For online instructions, access the link or scan this QR code:


<https://www.suport.ro/blog/instructiuni/masaj-picioare>


Package content

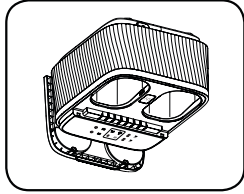1 x Foot massager

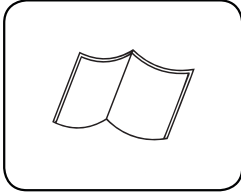1 x User manual

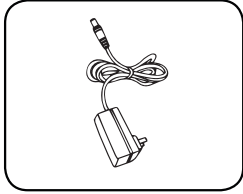1 x Power cord adapter

Thank you for buying this product!
We value your trust in Supporto.

Relax and improve your health with this
massage device!

“Take care of your body,
it's the only place you have to live.”

Jim Rohn

Device presentation

The aim of this product is to provide you with a relaxing and therapeutic
massage experience for your soles and feet, integrating advanced reflexology
techniques.

Benefits include deep muscle relaxation and reduction of accumulated
tension, improved blood circulation to the feet, improved immune system
functionality for disease prevention, maintenance of health and physical
well-being. UV sterilisation helps eliminate bacteria and germs for
increased hygiene, and heating adds a comforting and invigorating
sensation during use.

This machine is designed to bring massage, relaxation and revitalization right
into the comfort of your home.

Thank you	3
Device presentation	3
Importance of reading instructions	4
Package content	4
Ask your doctor advice before use	5
Safety rules and measures	5
General safety measures	5
Measures regarding the risk of electric shock	6
Fire risk instructions	6
Product description and parts	7
Correct use of the device	7
Setup	7
Operating instructions	8
Maintenance and storage	9
Troubleshoot	10
Technical specifications	11
Warranty terms	12
Recycling	12
Product certification	12

FOOT MASSAGER UV USER MANUAL

This product complies with European
conformity standards.

Made to order
Support Ortopedic SRL
djutor@supporto.ro
+40757246463
Importer: Support Ortopedic SRL