PISTOL DE MASAJ PROFESIONAL

MANUAL DE UTILIZARE

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
Appliance certification

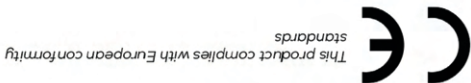Warranty terms

Legal provisions on consumer protection are applicable.

1. Warranty period for our products:

1.1. 2 years for individuals;

1.2. 1 year for legal entities.

2. The warranty is calculated as of the date of delivery of the product. During the warranty period, we will be responsible for repairing free of charge any defects caused by the quality of the product itself under normal conditions of use.

3. During the warranty period, if the defect is caused by the quality of the product, please contact us.

4. Please keep your purchase invoice, as this serves as proof of the product warranty.

5. Free repairs will not be carried out in the following cases:

5.1. Failure caused by not using the appliance in accordance with the instructions provided.

5.2. Damage caused by repairs by unauthorized persons.

5.3. Damage, scratches or breakage caused by moving or falling.

5.4. Damage caused by improper storage, maintenance or use by the user.

5.5. Damage caused by natural disasters.

~~Do not remove screws and do not try to disassemble the appliance, otherwise the product warranty becomes null.~~

Recycling

Please contact a local organization to recycle this product at the end of its useful life.
The product will be recycled after disposal of the batteries.

Precautions when disposing of batteries

The battery disposal procedure will be carried out by authorized personnel. The lithium-ion batteries used in this product are a valuable resource that can be recycled, so please dispose of them at the end of use by handing them in to a professional organization that can recycle rechargeable batteries. The device must be disconnected when the battery is removed. The battery must be disposed of safely. Please support the recycling of rechargeable batteries (lithium-ion batteries).

Battery removal procedure:

1. Use a professional tool to open the plastic housing on the non-key side.
2. Cut off the battery leads with scissors.
3. Remove the battery.

The battery used in this product is a special rechargeable battery and should not be used in other products.

- Do not overheat it by throwing it in the fire.
- Do not pierce with sharp objects, throw or make disassembly changes.
- Do not connect the positive and negative terminals of the battery directly to metal.
- Do not store batteries in the same place where you keep necklaces, hairpins.
- Do not charge, use or leave in hot places, such as near a fireplace or in the blazing sun, as these situations can cause hot spots, fires and breakdowns.
- Keep out-of-use rechargeable batteries away from babies and children to prevent accidental swallowing.
- In case of accidental swallowing, please seek medical attention immediately.
- In case of battery fluid leakage, do not touch the fluid with your hands or any other part of your body.
- If you come in to contact with battery fluid, wash immediately with plenty of water and seek emergency medical assistance.
- Failure to comply with the measures may result in fire, electric shock or personal injury.

Indications for use of massage heads

1. Spherical head: allows graduated intensity control, ideal for general use as well as sensitive areas: neck, elbow and knee area, back, forearm, arms, legs, shoulders.
2. Y-head – Suitable for large muscle groups: arms, legs, buttocks, calves, deep tissue.
- Warning! To be used along the spine only with professional guidance.
3. Pointed head – Used to precisely massage deep tissues, trigger points, paraspinal muscles, joint area, muscle knot release.
4. Round flat head – Used for medium and large muscle groups, muscle soreness, hamstring, biceps, calves.
5. Triangular head – Indicated for general use as well as for legs and buttocks.
6. T-head, mushroom – Suitable for flat muscles, thighs, abdomen.
7. Arched Head – Suitable for soft tissues, thighs, arms, calves, chest, abdomen.
8. Umbrella head with striations – Especially for anti-cellulite massage but also for legs, back, shoulders, soles and palms.

Maintenance and storage

Clean the device with a slightly damp cloth and wipe dry with a soft cloth. When storing or travelling, place it in its carrying case, away from liquids

Power, Battery and Charging

1. Power supply voltage: 5V/1A (slow charging) and 2A (fast charging).
2. Can be charged by phone charger (5V/1A or 5V/2A), external battery or laptop (slow).
3. Do not charge the device with very powerful chargers to avoid battery damage.
4. If it has not been charged for a long time, put it on charge to prolong battery life.
5. Do not keep in cold temperatures for a long time.
6. Charge the battery fully if it is to be stored for a long period.
7. Charge for six hours or more for a full charge before first use.
8. It is not recommended to fully discharge the battery.
9. Do not charge the device while it is in operation.
10. Do not leave the device charging unattended.

Technical specifications

Professional massage gun			Battery capacity	1500 mAh
Voltage			SV = 2A	Power rating
Message amplitude			6mm	Charging port
Charging time			2 hours	Duration of use
Gear speed			Smart Mode 1 - 99 levels	
Product size			150 mm * 58 mm * 189.5 mm	
Weight			675g	

Correct use of the device



Start up	Press the "⏻" button for 2 seconds to start the massage gun, initial speed is 1.
Speed setting	Once it is on, press "⊕", each press increases the speed level by one, the highest level is 99. Press "⊖" to decrease the speed.
Stop	Press "⏻" for 2 seconds to turn the device off.
Battery level display	Press "⊕" when switched off to display the current battery level.
Charging	When charging, the current charge level is displayed and when fully charged, the display shows 99.
Remember:	
1. The overlock(Turbo) function is activated when you put more pressure on a specific body part, its purpose is to raise the intensity of the massage quickly, whenever desired. The overlock function is turned off when the pressure is removed. 2. After 15 minutes of use, the massage gun will switch off automatically. If it is to be used further, it needs to be restarted.	

Gear light indicator

Product description and components

The massage gun is also known as a deep fascia relaxation device. It is a tool for repairing the soft tissues of the body. It has high torque, long range of action and adjustable massage speed. It is effective in reducing pain caused by lactic acid build-up following strenuous exercise and relieves discomfort in the muscles of people who have been sitting for a long time.

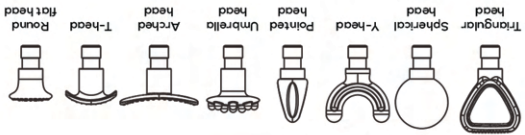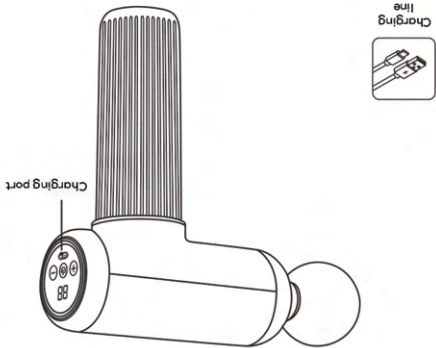Ask your doctor before use

Please consult your doctor if at least one of the situations below applies, but also in other situations if you are not sure whether the device is right for you.

1. Pregnant women, people suffering from diabetes and complications such as neuropathy, retinal damage.
2. Use of pacemaker, recent surgery, epilepsy, migraine, herniated disc, bone pushing forward.
3. Heart disease.
4. Recent joint replacement surgery or intravertebreal device placement, metal implants.

Safety rules and measures

To reduce the risk of electric shock, fire and personal injury, this product must be used in accordance with the following measures.

1. For adult use only.
2. Apply only light pressure when using, bruising may occur.
3. Use only on soft body tissues without causing pain or discomfort.
4. Interrupt use immediately if signs of pain or discomfort appear.
5. Carefully check that the equipment and battery level are normal before each use.
6. Keep your fingers, hair or other body parts away from the shaft and back of the vibrating head as there is a risk of pinching or catching them.
8. Do not immerse in water, keep it away from liquids or heat sources, and do not allow water to enter the ventilation holes or the appliance.
9. The device should not be used on the head or any hard or bony part of the body.
10. Use only on dry, clean skin. Apply gentle pressure to the skin.
11. Do not insert any objects into the ventilation holes of the device.
12. Do not clog the ventilation holes.
13. Move in each massage point for about 60 seconds.
14. Do not drop and do not use the device improperly.
15. The device is not intended for commercial use, it is only for individual use.
16. Do not leave packaging material within reach of children.
17. Do not use it on animals.
18. Do not repair the appliance yourself. Improper repair may expose the user to danger.
19. Do not use near fire or flammable substances.

Thank you for buying this product! We value your trust in Supporto.
Improve your health with this massage device!

**“Take care of your body,
it's the only place you have to live.”**
Jim Rohn

Contents

Contents	1
Importance of reading instructions	1
Thank you	2
Ask your doctor before use	3
Product description and components	3
Correct use of the device	4
Power, Battery and Charging	5
Technical specifications	6
Indications for use of massage heads	6
Maintenance and storage	7
Safety rules and measures	7
Recycling	8
Precautions when disposing of batteries	8
Appliance certification	9
Warranty terms	9

Importance of reading instructions

Taking the time to read these instructions is an investment in your safety and well-being. The massage gun is an electrical product! Read these instructions for use carefully! Failure to follow these instructions may result in property damage or personal injury. Keep the instructions for use and show them to other users. Hand this manual in when you dispose of the appliance.

For online instructions, click here or scan this QR code:

support.ro/blog/instruction/massage-gun

This product complies with European
conformity standards.

Made to order
Support Ortopedic SRL
Support Ortopedic SRL
+40757246463
ortopedic@supporto.ro
Importen Support Ortopedic SRL

USER MANUAL

PROFESSIONAL MASSAGE GUN

EN