RO

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## Warranty terms

Legal provisions on consumer protection are applicable.

1. Warranty period for this product:

1.1. 2 years for individuals.

1.2. 1 year for legal entities.

2. The warranty is calculated as of the date of delivery of the product.

During the warranty period, we will be responsible for repairing free of

charge any defects caused by the quality of the product itself under

normal conditions of use.

3. During the warranty period, if the defect is caused by the quality of the

product, please contact us.

4. Please keep your purchase invoice, as this serves as proof of the product

warranty.

5. The product loses its warranty in the following cases:

5.1. Failure caused by not using the appliance in accordance with the

instructions provided.

5.2. Damage caused by repairs by unauthorized persons.

5.3. Damage, scratches or breakage caused by moving or falling.

5.4. Damage caused by improper storage, maintenance or use by the user.

5.5. Damage caused by natural disasters or accidents.

5.6. Damage caused through the fault of the buyer (accidents including

transport)

5.7. Damage caused by falling, crushing, impact, exposure to water.

**Do not remove screws and do not try to disassemble the appliance, otherwise the product warranty becomes null!**

## Recycling

Don't throw the device away!

Please contact a local organization to recycle this product at the

end of its useful life.

## Product certification

This product complies with European conformity standards.

Product Code	REFONE	Message modes	3
Name	Shiatsu Foot Massager with heat	Mode of operation	Buttons
Voltage	12V  3000mA	Heating mode	Soothing heating
Power	36W	Timer	15 minutes
Voltage supply	100-240 V ~ 50/60 Hz	Dimension	34 x 35.5 x 7.6 cm
Weight	Aprox 1,93 kg	Adapter cable length	Aprox 1,8 m

## Technical specifications

Method of resolution	Reason	Error
<ol style="list-style-type: none"> <li>1. Plug the adapter into the socket.</li> <li>2. Check if the adapter plug is connected to the charging port of the device.</li> <li>3. Check that the socket has power.</li> </ol>	<p>Not connected to the socket</p>	<p>Device fails to start</p>
<p>Press the Power ON button (1)</p>	<p>Device in not ON</p>	<p>Device fails to start</p>
<ol style="list-style-type: none"> <li>1. If yes, allow the device to cool for 15 minutes before proceeding.</li> <li>2. The device may take longer to cool down when the ambient temperature is high.</li> </ol>	<p>Check if the device has been used for more than 30 minutes.</p>	<p>Device suddenly stops</p>
<p>Please contact your retailer or an authorised service centre if none of the above methods solve the problem.</p>	<p>Other reasons</p>	<p>There is a clattering noise during use</p>
<ol style="list-style-type: none"> <li>1. Continue using by pressing the "Power/ ON" button.</li> <li>2. After every 30 min of use switch off the appliance for 15 minutes.</li> </ol>	<p>The automatic switch-off function is activated after 15 minutes of operation.</p>	<p>This is normally caused by the operation of the motor and rollers during use.</p>

- Before cleaning, make sure the appliance is switched off and the adapter has been unplugged.
- Avoid immersing the appliance in any liquid during cleaning.
- Clean only with a soft, slightly damp sponge.
- Do not allow water or any other liquid to get inside the unit.
- Avoid using abrasive cleaners, brushes, gasoline, kerosene, glass/furniture polishes.

## Cleaning instructions

If you do not use the appliance for a long period of time, please keep it in a safe, dry and covered place to prevent dust accumulation, out of direct sunlight and out of reach of children.

Avoid contact with sharp edges or sharp objects that could cut or puncture the surface.

Please wrap and store the power cord under the bottom frame of the appliance. Do not hang the appliance from the power cord.

## Maintenance and storage

Turn off by pressing the power button on the control panel.  
Remove the plug from the power supply socket.

**The device will automatically switch off after 15 minutes (timer run out).**

 **Heating Only Mode:** Turn off the device, then press the Heating button to activate only the heating function.

Unit: Heating off  
Dim red: Gentle heating  
Bright red: Strong heating  
Press to activate heating. (default high level).

### Heat Button

Red - High intensity  
Yellow - Medium intensity  
Green - Low intensity

### Intensity Button

Press to change intensity:


## Instructions for use

Choose a comfortable place to sit (armchair, sofa, chair) and place your feet on the massage roller surface. Make sure you are sitting in a comfortable position, without any tension.

Do not stand on the device, as this will break it!

Avoid using the device in bathrooms or damp areas.

Use the panel below to customize your massage experience.

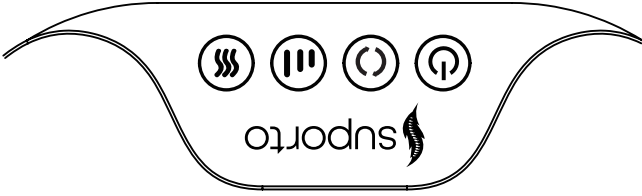**Power Button**

Press to turn on the device, it starts to work in Mode 1, press again to turn off the device

**Mode Button**

**Mode 1 – Automatic (Yellow LED)**

The upper and lower nodes move to the right, and the middle node moves to the left, then reverses. Variable intensity.


**Mode 2 – Green LED**

The top and bottom nodes move to the right, and the middle node moves to the left.

**Mode 3 – Red LED**

The top and bottom nodes move to the left, and the middle node moves to the right.

## Product description and parts
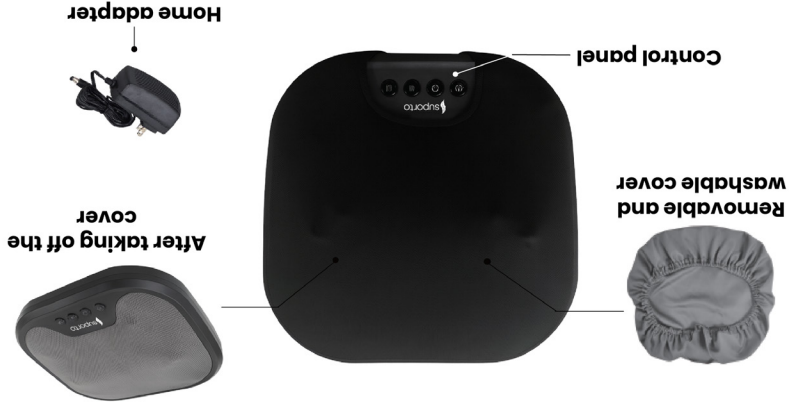
## Correct use of the device

Make sure all the components are free from visible defects before use.

If in doubt, do not use the device and contact the retailer or an authorized

service center.

Follow the instructions described below only after you have read, understood and observed all the safety rules and measures described in this document.

### Installation

1. Place the device where you want to use it. The surface must be flat and dry.
2. Connect the cable from the adapter to the power port on the device.
3. Plug the adapter into an electrical outlet (the adapter is safe for use with 110-120V & 220-240V).
4. After plugging it in, operate the controls using the buttons as described below.

9. Stop using the product if you notice malfunction or any other abnormal operation.

### **Measures regarding the risk of electric shock**

The appliance is very well insulated against the risk of electric shock.

Since it is an electrical device, it shall be used with care and caution to avoid the risk of electric shock:

1. Do NOT use in the following cases:

- 1.1. At a mains voltage different from that stated on the appliance.
- 1.2. If the appliance or its accessories have any noticeable defects (burningsmell, broken or damaged wire, visible burns, physical defects, etc.).
- 1.3. During storms.
- 1.4. In the presence of sources of increased humidity (bath, sauna).
- 1.5. If water has entered the appliance, the power cord, power plug or anywhere inside the product.

2. Do not pull the power cord or the appliance to unplug the plug. Never hold the appliance by the power cord.
3. Keep the power cords away from hot surfaces.

4. Handle only with dry hands. Do not unplug the power plug with wet hands. Do not use the appliance for more than 30 minutes consecutively, there is a danger of overheating. Allow the appliance to cool for at least 30 minutes after every hour of use.
6. If the appliance has been dropped or otherwise damaged as a result of being dropped, it may be used only after being checked by authorized personnel.
7. Do not repair the appliance yourself, contact the retailer or an authorized service center. Improper repairs may expose the user to serious risk of injury and the product warranty can become null.

8. Do not pull, twist or bend the power cord. Avoid placing the power cord over or in contact with sharp or pointed objects.

### **Fire risk instructions**

Do not use in the following cases:

1. Near fire or a heat source.
2. Near flammable or corrosive substances.
3. During sleep. Be careful not to fall asleep while the machine is still working.

## Ask your doctor advice before use

1. Ask your doctor advice before use if you are in at least one of the following cases:  
People with chronic foot conditions, who have recently undergone surgery on the feet, or who have open wounds, swelling, trauma, or burns to the feet or lower limbs.
  2. People with diabetes.
  3. People with pain whose cause is unknown.
  4. People with psychiatric conditions.
  5. People with pacemakers or implanted artificial organs.
  6. People with heart dysfunction, tumors and hypertension.
  7. People with brain haemorrhage and thrombosis.
  8. People with haemorrhagic diseases, septicæmia, fever, etc.
  9. Women who are pregnant or menstruating.
  10. People with advanced osteoporosis or fractures of the legs and back.
  11. People with heat intolerance.
- Please also consult your doctor in other cases if you are not sure whether the device is right for you!
- Avoid using the product before bedtime, it can affect your sleep.
- If you experience pain or discomfort during the session, stop using the device and consult your doctor!

## Safety rules and measures

### General safety measures

1. This appliance is NOT a toy. Close supervision is required when this appliance is used by or near children or persons with disabilities. Do not allow children to play with the appliance.
2. To avoid excessive stimulation of muscles and nerves, the recommended duration of massage should not exceed 15 minutes per session. Allow a waiting time between sessions.
3. Do not use the machine in damp or dusty environments as this may lead to malfunction or electric shock.
4. The appliance is not intended for commercial use, but for individual use only.
5. Do not use the appliance in the bathroom or on the bed to avoid risks. Be sure to store it in a safe place away from the edge of bathtubs or sinks where it could fall.
6. Do not sit on the appliance during use.
7. Protect from excessive heat or direct sunlight to prevent damage.
8. Do not use the appliance if the fabric cover is damaged.


## Importance of reading instructions

Taking the time to read these instructions is an investment in your safety and well-being.

The back massager is an electrical product! Read these operating instructions carefully!


Failure to observe the following instructions may result in property damage and personal injury. Keep the operating instructions and present them to other users. Hand in this manual when you dispose of the appliance.

For instructions in online version, go to this link or scan the QR code:

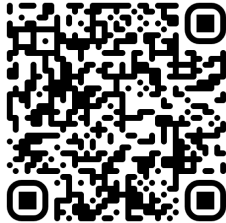
<https://www.supporto.ro/blog/instructiuni/masaj-picloare>


## Package content

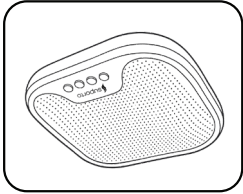1 x Back massager + removable cover

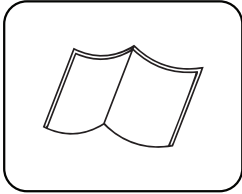1 x User manual

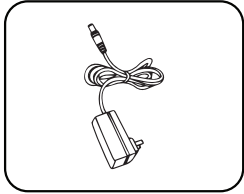1 x Power cord adapter

This foot massager helps you unwind and recharge after a long day with its multiple massage functions and adjustable intensity levels. The soothing infrared heat boosts circulation and eases muscle tension, while the built-in 15-minute auto-timer ensures comfortable use.

Compact, easy to use, and simple to maintain, it's the perfect way to enjoy moments of comfort and relaxation right at home.

## Device presentation

“Take care of your body,  
it's the only place you have to live.”

Jim Rohn

Thank you for buying this product!  
We value your trust in Supporto.  
Relax and improve your health with this  
massage device!

# Contents

Thank you .....	3
Device presentation.....	3
Importance of reading instructions.....	4
Package content.....	4
Ask your doctor advice before use .....	5
Safety rules and measures.....	5
General safety measures .....	5
Measures regarding the risk of electric shock.....	6
Fire risk instructions .....	6
Correct use of the device .....	7
Setup.....	7
Operating instructions.....	8
Maintenance and storage.....	10
Troubleshoot.....	11
Technical specifications.....	12
Warranty terms .....	13
Recycling.....	13
Product certification.....	13

# REFLEXO ONE - SHIATSU FOOT MASSAGER WITH HEAT

Reflexo One - Shiatsu Foot  
Massager with heat

This product complies with European  
conformity standards.

Made to order  
Support Ortopedic SRL  
Support Ortopedic SRL  
+40757246463  
supportor@supto.ro  
Importeur: Support Ortopedic SRL  
www.supto.ro

supporto 

EN