


LASCOD

If you like Kromopan you will love KromopanSil

New generation of
Addition Curing Silicones

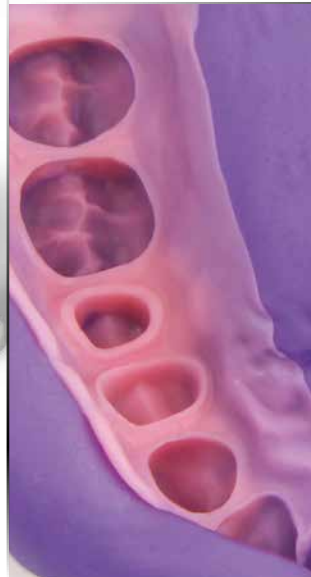

KromopanSil


addition curing silicone

Polyvinylsiloxane addition curing silicone is a versatile and efficient product guaranteeing maximum precision in humid conditions. KromopanSil has different grades of hardness (Hard/Soft) and consistency (Superlight/Light/Regular/Heavy) which attain the best results adapting to every impression technique.



A long working time increases the product workability. At the same time a short setting time in mouth is guarantee for the best patient's comfort.



KromopanSil is highly thixotropic and hydrophilic. The material stays in place after dispensing and its flows well under pressure. KromopanSil easily adapts to the anatomical shape of dental structures and soft tissues, thus reproducing perfectly integral margins.

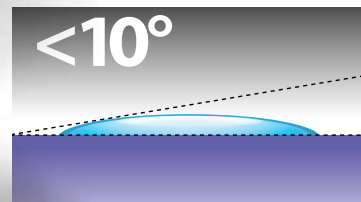

The contact angle measurement (wettability index) is less than 10°.

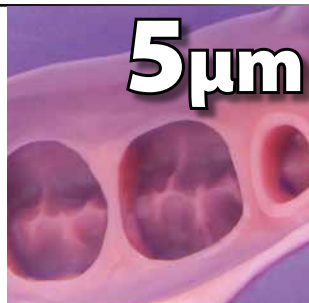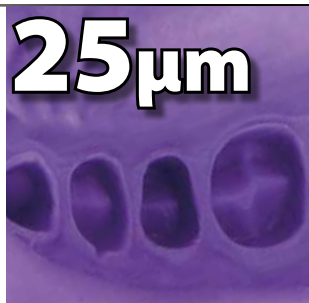

	MIXING TIP COLOURS		
SUPERLIGHT	Yellow		
LIGHT	Yellow		
REGULAR		Pink	
HEAVY			Teal
MONO		Pink	
BITE			Teal



KromopanSil cartridges perfectly match our mixing tips. Mixing tips are available in 48pcs packs, easily recognizable by colours: yellow (GNS200), pink (GNS250) and teal (GNS270).

KromopanSil is available in two different setting releases: normal and fast version, the fast set is faster both in the working time and the setting time in mouth.

	Working time including mixing time (23°)		Minimum time in mouth (37°)	
	NORMAL SET	FAST SET	NORMAL SET	FAST SET
PUTTY SOFT	2' 30"	1' 15"	2' 15"	1' 30"
PUTTY HARD	2' 30"	1' 15"	2' 15"	1' 30"
SUPERLIGHT	2'	1' 20"	1' 30"	60"
LIGHT	2'	1' 20"	1' 30"	60"
REGULAR	2'	1' 20"	1' 30"	60"
HEAVY	2'	1' 20"	1' 30"	60"
MONO	2'	-	1' 30"	-
BITE	-	45"	-	60"

KromopanSil is also available in kit packs containing 2x150ml putty base and catalyst, 1x50ml cartridge, 6 mixing tips and 6 oral tips.



KromopanSil kit	NORMAL	FAST
PUTTY HARD + LIGHT BODY	KRS080	KRS580
PUTTY HARD + SUPERLIGHT BODY	KRS081	KRS581
PUTTY SOFT + LIGHT BODY	KRS082	KRS582
PUTTY SOFT + SUPERLIGHT BODY	KRS083	KRS583

KromopanSil


MONO

monophase addition curing silicone

KromopanSil Mono, developed by Lascod, is an impression material with a medium viscosity, ideal for monophase technique. Recommended for crowns and bridges impression taking, implantoprosthesis, inlay, onlay, partial and total edentulism. Thanks to its characteristic viscosity it can be used directly on the surfaces to be reproduced and on impression tray. It's syringable through the included mixing tips or through the specific syringe for elastomers. Given its optimal flexibility **KromopanSil Mono** can be easily removed by the mouth without structural alteration.

KromopanSil Mono can be used also with Super Light Body / Light Body / Regular Body wash viscosity as support material inside individual impression tray or standard stainless steel impression tray.

- Ideal viscosity for monophase technique
- High details definition (5 µm)
- Maximum thixotropicity and hydrophilicity
- Fast setting time in mouth
- High resistance to tear
- High elastic recovery
- High dimensional stability over time
- Slight mint flavour
- Resistant to aseptic cleaners

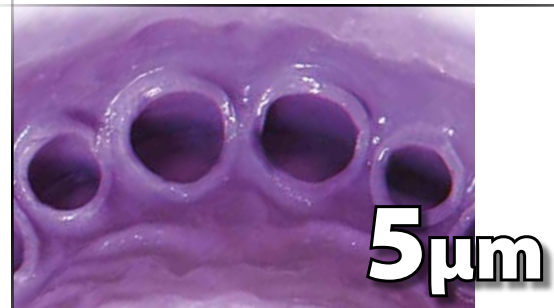

KromopanSil

silicone for bite registration

BITE

KromopanSil Bite is an addition silicone (polyvinylsiloxane) made for orthodontic bite registration, registration keys for gnathologic registrations, inter-maxillary registration keys for centered positions, eruptions and ectopic eruptions, registration for cephalometric analysis, which requires a CAD scanning (optical, laser and tactile).

KromopanSil Bite characteristics grant a working base which allow the technicians to obtain the maximum precision.

The high dimensional stability over time ensures optimal conservation even in the presence of temperature changes.












The useful hardness degree (**95 Shore A**) achieved in short time, allows to work on masses widely stable facilitating the finishing stages by cutters or burs.

- Easily and safely workable thanks to double cartridge system 1:1
- Maximum control during positioning
- Thixotropic
- Imperceptible viscosity for patient
- Exact reproduction of occlusal details (5µm)
- Short setting time
- Minimum dimensional variation
- Undeformable / stable over time
- Easy insertion and removal from the model
- Easy finishing by burs

packaging (FOR KITS LOOK AT THE PREVIOUS PAGE)

*** ALSO AVAILABLE IN ECONOMY PACK**
4 x 50 ml cartridges
without mixing tips

	KromopanSil	NORMAL	FAST	
	PUTTY HARD	KRS020	KRS520	2 x 300 ml base + catalyst
	PUTTY SOFT	KRS010	KRS510	2 x 300 ml base + catalyst
	HEAVY BODY*	KRS070	KRS570	2 x 50 ml cartridges + 12 mixing tips
	SUPERLIGHT BODY*	KRS040	KRS540	2 x 50 ml cartridges + 12 mixing tips
	LIGHT BODY*	KRS050	KRS550	2 x 50 ml cartridges + 12 mixing tips
	REGULAR BODY*	KRS060	KRS560	2 x 50 ml cartridges + 12 mixing tips
	MONO* (50 ml)	KRS400	-	2 x 50 ml cartridges + 6 mixing tips
	MONO (380 ml)	KRS480	-	2 x 380 ml cartridges + 2 mixing tips
	KromopanSil BITE			
	<i>Silicone for bite registration *</i>	-	KRS095	2 x 50 ml cartridges + 12 mixing tips

The photographic images here reproduced are purely indicative and are not necessarily identical to the actual products.



LASCOD S.p.A. - Via L. Longo, 18 - 50019 Sesto Fiorentino, Florence, Italy - tel.: +39 055 4215768 - lascod.italy@lascod.it
www.lascod.com

FOLLOW US!



KRSGB4 - 02/2020